TERMOMETRI DIGITALI

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Hangzhou Sejoy Electronics & Instruments Co. Ltd
Building 2, No.202, Zhenzhong Road, West Lake Economy & Technology Zone 310030 Hangzhou China

Shanghai International Holding Corp. GmbH (Hamburg) - Eiffelstrasse 80, 20537 Hamburg, Germany

Distrib. MORETTI S.P.A. Via Bruxelles, 3 - Melegnano 52022 Cavriola (Arezzo) - Italy - www.morettispa.com MADE IN P.R.C.

MA RM305_310_320_01A_03-2016

DIGITAL THERMOMETER

- IIa Class Medical Device
According to Directive 93/42/EEC and further modifications

ENGLISH

1. CODES

- RM305 Rigid digital thermometer, 60 sec - Water proof
- RM310 Flex digital thermometer, 60 sec - Water proof
- RM320 Flex digital thermometer, 30 sec - Water proof - Display Jumbo

2. INTRODUCTION

Thank you for purchasing GRADO a by Moretti Thermometer. This user's manual provides some suggestions as how to correctly use the product you have chosen and gives some valuable advice for your safety. Please read through the manual carefully before using the product. Should you have any queries, please contact your retailer for advice and assistance. In case of doubts, we recommend to contact the dealer who can help and advice you properly.

3. INTENDED USE

The thermometer is a device designed to measure body temperature via oral, axillary or rectal.

CAUTION!

- Do not use the product for a purpose not indicated in this manual.
- Moretti S.p.A declines all responsibilities for any consequences resulting from an incorrect use of this product and from unauthorized alteration to the frame of the product.
- The manufacturer reserves the right to change the information contained in this document without previous notice

4. DECLARATION OF CONFORMITY

Moretti S.p.A. declares under his sole responsibility that the product placed in the market by the same Moretti S.p.A. is in conformity to the relevant provisions of Directive 93/42/CEE on medical devices and subsequent amendments and are produced by Hangzhou Sejoy Electronics & Instruments Co. Ltd according to the standards indicated in the regulations.

- CEI 60601-1 - Part 1: General requirements for basic safety and essential performance.
- EN 60601-1-2 - Part 1: General requirements for safety - Collateral standard: Electromagnetic compatibility.
- ISO 80601-2-56:2009: Medical electrical device - (Part 2-56: Particular requirements for basic safety and essential performance of clinical thermometers for measuring body temperature).

5. GENERAL WARNINGS

- Read and understanding this instruction manual.
- Choking Hazard: Thermometer cap and battery may be fatal if swallowed. Do not allow children to use this device without parental supervision.
- Do not store or keep in direct sunlight, or in high temperature.
- Shelf life - the thermometer use limit is defined by the wear of the parts.
- Do not use in the ear. Designed use is for oral, rectal, and armpit (axilla) readings only.
- Do not place thermometer battery near extreme heat as it may explode.
- Do not allow children to play on or operate the thermometer.
- It is recommended the performance should be checked every two years.
- Use of the probe cover may result in a 0.1°C discrepancy from actual temperature.
- Remove battery from the device when not in operation for a long time.
- The use of temperature readings for self-diagnosis may lead to the worsening of existing disease conditions.
- Do not attempt measurements when the thermometer is wet as inaccurate readings may result.
- Do not bite the thermometer. Doing so may lead to breakage and/or injury.
- Do not attempt to disassemble or repair the thermometer. Doing so may result in inaccurate readings.
- After each use, disinfect the thermometer especially in case the device is used by more than one person.
- Do not force the thermometer into the rectum. Stop insert and abort the measurement when pain is present. Failure to do so may lead to injury.
- Do not use thermometer orally after being used rectally
- For children who are two years old or younger, please do not use the devices orally.

6. ELECTROMAGNETIC COMPATIBILITY WARNINGS

This equipment has been tested and it complies with medical devices limits belonging to EN 60601-1-2:2007. These limits have been designed to provide reasonable protection against harmful interference in a typical medical installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to other devices in the vicinity. However, there is no guaranteed that the interference will not occur in a particular installation. If this equipment does cause harmful interference to other devices, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference following one or more of the suggestions below:


- Reorient or relocate the receiving device.
- Increase the separation between the equipment.
- Consult the manufacturer or field service technician for help.

CAUTION!

Do not use in proximity of the device, mobile phones and other devices, which can generate strong electrical or electromagnetic fields as they can interfere with the operation of the device.

7. SIMBOLOGY

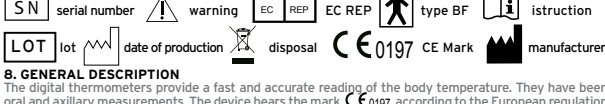8. GENERAL DESCRIPTION

The digital thermometers provide a fast and accurate reading of the body temperature. They have been designed for rectal, oral and axillary measurements. The device bears the mark CE 0197 according to the European Regulation for Medical Devices. The device isn't subject to the calibration obligation though a quality control in accordance with directive 93/42/CEE.

9. WARNINGS FOR USE

CAUTION!

- The device performance can be compromised if one or more of the following events occur:
 - Operation outside the manufacturer's stated temperature and humidity range
 - Storage outside the manufacturer's stated temperature and humidity range
 - Patient temperature is below ambient temperature

- Do not bend, drop the thermometer.
- Do not boil probe to disinfect.
- Do not disassemble, except in case of battery replacement.
- The thermometer should only be used under the supervision of an adult. Prohibit children from biting, hitting and removing the battery cover; dislodging and potentially swallowing the battery.
- Do not use the device in case of strong electromagnetic field such as microwave oven or cell phone.

Body temperature

- Body temperature can vary from 36.1°C to 37.2°C. The generally accepted "normal" temperature is 37.0°C.
- Usually the temperature is lower in the morning than in the afternoon. The temperature "normal" varies depending on the environment where the measurement is made. While the oral temperature follows the above written, the rectal temperature is generally higher than 0.5 °C; on the contrary the axillary temperature is less than 0.5 °C.

Note: Temperature readings are available in the Fahrenheit or Celsius scale (°C/°F); located in the upper right corner of LCD.) With the unit off, press and hold the On/Off button for approximately 4 seconds to change the current setting.

10. HOW TO USE

Measure the temperature

- Press On/Off button.
- You will hear a beep and the screen will display "188,8E", then will display the last reading and then the auto-test temperature.
- Finally the screen will display "Lo °C" and the symbol "°C" start flashing.
- The thermometer will start to measure the temperature as soon as it is exposed to body heat.
- The °C symbol will continue to flash until completion of the temperature taking process.
- After the measurement, the thermometer will emit a sequence of acoustic signals.
- Push On/Off button to turn the unit off. If you do not turn off the unit, it will switch off automatically.
- Store the thermometer in safety place, as manual is written.

Note: If the measured temperature is less than 32°C the screen will display Lo, if the measured temperature is higher than 42.9°C the screen will display Hi. Alarm beeps more rapidly when temperature reaches 37.8°C (100°F) or higher.

Method of measurement

There are three methods of measuring body temperature. Slight temperature differences may occur depending on the method of measurement.

- Oral measurement:** using the diagram, place the probe tip under the tongue as near as possible to the heat pocket. Close your mouth and breathe evenly through the nose to prevent the measurement from being influenced. Normal temperature between i 35.7°C and 37.3°C.
- Axillary measurement:** the thermometer should be placed under the arm so the probe tip is touching the skin. Hold the thermometer parallel to the body and position the arm on the chest so the tip is well covered by the arm. This insures the probe is not affected by the room's air. Normal temperature between i 35.2°C and 36.7°C.
- Rectal measurement:** rectal temperature should be taken when is impossible or impractical to take oral or axillary temperature. The probe should be inserted a maximum of 1-2 cm. Normal temperature between i 36.2°C and 37.7°C.

Replace the battery

- The battery should be replaced when the screen display Shows "Bat".
- Pull battery cover off as shown in Fig. below and use a screwdriver.
- Use pointed object such as a pen to remove old battery. Replace with new 1.5DC button type.
- Position the battery cover on the thermometer.

illustrative purposes only

11. MAINTENANCE

The GRADO devices by Moretti S.p.A are checked carefully and supplied with EC mark, once launched on the market. For patient and physician's safety we recommend to let check the device by the manufacturer or by laboratory authorized every 2 years. In case of repair, use only original spare parts and accessories.

12. CLEANING AND DISINFECTION

Cleaning

Wipe the thermometer with a soft clean cloth. For stubborn stains, wipe the thermometer with a cloth that has been dampened with water or a neutral detergent solution and then wring thoroughly. Finish by wiping with a soft dry cloth.

Observe the following to prevent damage to the thermometer:

- Do not use benzene, thinner, gasoline or other strong solvents to clean the thermometer
- Do not soak the sensing section in alcohol for long periods of time or attempt to sterilize it using hot water at 50°C or higher
- Do not use ultrasonic washing to clean the thermometer

Disinfection

If you need to disinfect the device use isopropyl rubbing alcohol.

Note: Never use acids, alkalis or solvents such as acetone.

13. CONDITIONS OF DISPOSAL

General conditions of disposal

In case of disposal don't use the inserting container for municipal waste. We recommend to dispose the thermometer in the appropriate disposal areas for recycling.

Information for the correct disposal of the product in accordance with European Directive 2002/96/CE:

At the end of its life, the product must not be disposed of along with other domestic waste. The users must dispose of this equipment by bringing it to a specific recycling point for electric and electronic equipment or at retailers that provide this service. By ensuring these batteries are disposed of correctly, you will help prevent potentially negative consequences for the environment and human health which could otherwise be caused by inappropriate waste handling of the battery. To remark the need to dispose of electrical equipment separately, the products is marked with crossed mobile waste bin the product bears the barred bin symbol.

Disposal of waste batteries (directive 2006/66/EC):

This symbol on the battery or on the packaging indicates that the battery provided with this product shall not be treated as household waste. By ensuring these batteries are disposed of correctly, you will help prevent potentially negative consequences for the environment and human health which could otherwise be caused by inappropriate waste handling of the battery. The recycling of the materials will help to conserve natural resources. At the end of their life hand the batteries over to the applicable collection points for the recycling of waste batteries. For more detailed information about recycling of this product or battery, please contact your local Civic Office, your household waste disposal service or the shop where you purchased the product.

14. SPARE PARTS AND ACCESSORIES

For spare parts and accessories refer to main catalogue.

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EMISSION TEST	COMPLIANCE	ELECTROMAGNETIC ENVIRONMENT-GUIDANCE
RF emissions CISPR 11	Group 1	The device uses RF Energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment
RF emissions CISPR11	Class B	The device is suitable for use in all establishment. Including domestic establishments and those directly connected to the public low-voltage power supply network
Harmonic emissions IEC61000-3-2	N/A	
Voltage fluctuations/flicker emissions IEC61000-3-3	N/A	

Guidance and manufacturer's declaration - Electromagnetic immunity -

The devices are intended for use in the electromagnetic environment specified below. The user of this device should make sure it is used in such an environment.

Immunity Test	IEC60601 test level	Compliance	Electromagnetic Environment-Guidance
Electrostatic Discharge (ESD) IEC61000-4-2	± 6kV contact ± 8kV air	± 6kV contact ± 8kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/ burst IEC61000-4-4	± 2kV for power supply line ± 1kV for input/output line	N/A	
Surge IEC61000-4-5	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	N/A	
Voltage dips, short interruptions and voltage variations on power supply input lines IEC61000-4-11	<5 % UT (>95 % dip in UT)for 0,5 cycle 40 % UT (60 % dip in UT)for 5 cycles 70 % UT (30 % dip in UT)for 25 cycles <5 % UT (>95 % dip in UT)for 5 sec	N/A	
Power frequency (50/60Hz) magnetic field IEC61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of atypical location in a typical commercial or hospital environment.

Note: UT is the a.c. mains voltage prior to the application of the test level

Guidance and Manufacturer's Declaration - Electromagnetic Immunity -

The devices are intended for use in the electromagnetic environment specified below. The user of this device should make sure it is used in such an environment.

Immunity Test	IEC60601 test level	Compliance	Electromagnetic Environment-Guidance
Conducted RF IEC 61000-4-6	3Vrms150 kHz to 80 MHz outside ISM banda	N/A	Portable and mobile RF communications equipment should be used no closer to any part of this device, including cables, than there commended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance d = 1.2√ P d = 1.2√ P 80 MHz to 800 MHz d = 2.3√ P 80 MHz to 2.5G MHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).b Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey c, should be less than the compliance level in each frequency ranged. Interference may occur in the vicinity of equipment marked with the following symbol:
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

NOTE 3 Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and landmobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the device

NOTE 4 Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 KHz to a 80MHz d=1.2√ P	80MHz to 800MHz d=1.2√ P	800MHz to 2.5 GHz d=2.3√ P
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

The user can prevent electromagnetic interference by maintaining a minimum distance between the thermometer and the device that might causes interference.

16. TECHNICAL FEATURES

- WEIGHT (INCLUDING BATTERY): RM305 10gr; RM310 10gr; RM320 18gr
- RANGE OF TEMPERATURE: 32.0°C ± 0.1°C - 42.9°C ± 0.1°C
- ACCURACY: ± 0.2°C below i 35.5°C; ± 0.2°C over i 42°C
- MEMORY: For storing the last measured value
- BATTERY: One 1.5V DC button battery
- BATTERY LIFE: Approx 200 hours
- DISPLAY: Liquid crystal display (LCD) 3 1/2 DIGITS
- USING ENVIRONMENT: 5°C to 40°C/ 15% to 95% RH
- STORAGE AND TRANSPORTATION CONDITION: -20°C to 55°C/ 15% to 95% RH
- CLASSIFICATION: Type BF
- DIMENSION (LxWxH): RM305 (13,1x1,9x1,0cm); RM310 (13,9x2,7x1,2cm); RM320 (13,4x3,1x1,65cm)
- INGRESS PROTECTION RATING: RM305 IP27; RM310 IP27; RM320 IP27
- REGULATION: EN 12470-3; ISO 80601-2-56:2009

17. WARRANTY

Every Moretti product is guaranteed for 2 (two) YEARS from the purchase date, except the restrictions specified below. The warranty shall not apply in case of damage caused by improper use, abuse, alterations of the product or any other use not in compliance with the user manual.

The proper use of the assembly is indicated in the user manual. Moretti S.p.A. shall not be liable for any damages or personal injuries due to incorrect assembly/improper use of the device, not compliant with the user manual. Moretti does not guarantee its products against damages or faults in the following cases: natural disasters, non authorized maintenance or repairs, damages caused by electrical supply problems (if provided), use of parts or components not provided by Moretti, in accordance with the guidelines and instructions for use, non authorized alterations, damages due to the transport (different from Moretti original delivery) or failure in maintaining the device as indicated in the manual. The components subject to wear are not covered by this warranty if the damage is due to the normal use of the product

Rechargeable batteries warranty (If provided)

The original batteries and the spare batteries are covered by a 90 days warranty during a period of 6 months, with regards to manufacturing defect or to what concern the law regulation. If charge batteries are unused for more than 3 consecutive months, the warranty is void. If exhausted batteries are unused for more than 3 consecutive days, the warranty is void.

18. REPAIRING

Warranty repair

If a Moretti product shows material or manufacturing faults during the warranty period, Moretti S.p.A., together with the client, will consider whether the product fault is covered by warranty. Moretti S.p.A., at its unquestionable discretion, may replace or repair the item under warranty in a Moretti reseller point or at its own qualified offices. Labour costs for repairing the product may be at Moretti's expense if it is determined the reparation is under warranty. A reparation or replacement does not renew nor extend the warranty.

Repair not covered by warranty

A product, not covered by warranty, may be returned for reparation only after receiving prior authorization from Moretti client service. Labour and delivery costs for an out-of-warranty reparation will be fully at the client's or reseller's expense. Reparations of out-of-warranty products are guaranteed for 6 (six) months, from the day the client receives the repaired product.

Non defective device

The client is informed if, after examining and testing a returning product, Moretti decides the product is not defective. The product will be sent back to the client and all the return delivery costs will be at his expense.

19. SPARE PARTS

The original Moretti spare parts are guaranteed for 6 (six) months from the day you receive them.

20. EXEMPT CLAUSES

Except what clearly specified in this warranty and in accordance with the law, Moretti S.p.A. does not offer its other declaration, warranty or condition, expressed or implicit, including any declaration, warranty or condition of merchantability, of non violation or non interference, of suitability for a specific aim. Moretti S.p.A. does not guarantee that the use of Moretti products will be continuous and faultless. The duration of possible implicit warranties that may be applied in accordance with the law is limited to the warranty period.

Some states or countries do not allow limitations to the duration of an implicit warranty, or the exoneration, or the limitations of accidental or indirect damages related to products for consumers. In such states or countries some exemptions or limitation of this warranty may not be applied to the user. This warranty is subject to modifications without any prior notice.

Hangzhou Sejoy Electronics & Instruments Co. Ltd
Building 2, No.202, Zhenzhong Road, West Lake Science & Technology Zone 310030 Hangzhou China

Shanghai International Holding Corp. GmbH (Hamburg) - Eiffestrasse 80, 20537 Hamburg, Germany

Distrib. MORETTI S.P.A. Via Bruxelles, 3 - Melegnano 20122 Cavriaglia (Arezzo) - Italy - www.morettispa.com MADE IN P.R.C.

MA RM305_310_320 01 A_03-2016