**04550** rev.001  
HOME WALL BARS  
SPALLIERA CASA

---

# HOME WALL BARS

## SPALLIERA CASA

### Use and maintenance manual

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**CHINESPORT SPA**  
Via Croazia, 2  
33100 – Udine – IT  
Tel: +39 0432 621 621  
Fax: +39 0432 621 620  
[www.chinesport.it](http://www.chinesport.it)  
[info@chinesport.it](mailto:info@chinesport.it)

ENGLISH .....	3
1 GENERAL INFORMATIONS .....	3
1.1 SYMBOLS .....	3
2 ASSEMBLY AND INSTALLATION .....	4
2.1 LIST OF PARTS .....	5
3 ASSEMBLY.....	5
4 MAINTENCE.....	6
4.1 ROUTINE MAINTENANCE .....	6
5 CLEANING AND DISINFECTION .....	7
5.1 WASHING.....	7
5.2 DISINFECTION.....	7
6 WARRANTY .....	7
7 MAINTENANCE RECORD SHEET .....	8
ITALIANO .....	9
1 INFORMAZIONI GENERALI .....	9
1.1 SIMBOLI .....	10
2 MONTAGGIO ED INSTALLAZIONE .....	10
2.1 CONTENUTO DELLA CONFEZIONE .....	11
3 ASSEMBLY.....	11
4 MANUTENZIONE .....	13
4.1 MANUTENZIONE ORDINARIA.....	13
4.2 MANUTENZIONE STRAORDINARIA .....	13
5 PULIZIA E DISINFEZIONE.....	13
5.1 LAVAGGIO .....	13
5.2 DISINFEZIONE.....	13
6 GARANZIA.....	13
7 REGISTRAZIONE INTERVENTI.....	14

## 1 GENERAL INFORMATIONS

<b>PURPOSE AND CONTENT OF INSTRUCTION</b>	<p>The purpose of this manual is to provide users with the necessary information to allow proper use of the equipment and make it possible to use it independently and safely. This manual contains information on technical, operational and maintenance aspects, as well as spare parts and safety. It is also intended for technicians responsible for maintaining the equipment. This Instruction Manual must always be kept with the equipment and must be carefully preserved in order to avoid compromising its readability. If the equipment is sold to third parties this manual has to be handed over with the other documents as required by European Directive 93/42 EEC.</p> <p><u>This manual is not a medical handbook</u></p>
<b>MANUFACTURER</b>	<p><b>CHINESPORT S.P.A</b>, Via Croazia, 2-33100 Udine, Italia. Phone. +39 0432 621 621 – Fax+39 0432 621 620 – web: <a href="http://www.chinesport.com">www.chinesport.com</a> The company produces according to UNI EN ISO 13485:2012.</p>
<b>INTENDED USE AND EQUIPMENT DESCRIPTION</b>	<p>Swedish beech-wood wall bars with single span, are natural varnished and supplied with anchoring white steel clamps galvanized. The white steel clamps have to be used at the four angles to fix the wall bars. The wall bars are made up of 12 cylindrical rungs. They have a diameter of 3,3 cm and there is a space of 14 cm between each other</p>
<b>FURTHER INFORMATIONS AND DOCUMENTATION</b>	<p>Please visit <a href="http://www.chinesport.com">www.chinesport.com</a></p>
<b>STORAGE</b>	<p>In case of storage the following ambient conditions must be respected: relative umidity 10% / 90 %; temperature -10°C / +50°C.</p>
<b>SIDE EFFECTS AND CONTRAINDICATIONS</b>	<p>There is no collateral effects or contraindications during proper use. Use the device always in the proper way. The knowledge of this instruction manual is recommended.</p>
 <b>SPECIAL WARNINGS AND USE RESTRICTIONS</b>	<p>These devices must not be used:</p> <ul style="list-style-type: none"><li>• without making all the settings required and recommended setups.</li><li>• when use is not advisable or causes the patient serious discomfort.</li><li>• if the original product has been tampered with and/or modified</li><li>• when setted incorrectly or when there is any suspicion that it may harm the patient</li><li>• with patients over 120 kg</li><li>• with more than 1 patient</li></ul> <p>The plastic and wood parts can catch fire if brought into contact with naked flames. Do not spill liquids on the metal parts. The physical and psychological suitability of the user to make use of the aids must be approved by a specialist.</p>

### 1.1 SYMBOLS

The symbols found in this manual and on the machine are used to draw the user's attention to important safety warnings and useful information. Observe all warning, caution and danger messages.

	<b>WARNING</b> This warning symbol identifies operating/maintenance instructions or procedures which, if not observed, could result in injury or danger to the user.
	<b>ELECTRIC SHOCK HAZARD</b> This warning symbol identifies operating/maintenance instructions or procedures which, if not observed, could result in a risk of electric shock.
	<b>CAUTION</b> A Caution symbol is used to warn the user that failure to observe a certain operating or maintenance procedure, general rule or condition, could result in damage to the equipment.
	<b>NOTE</b> Notes draw the user's attention to key information on topics, applications or procedures.
	<b>TIPS</b> Tips draw the user's attention on the quickest or easiest way to do something.
	<b>WORKERS</b> This symbol indicates the minimum number of people recommended for installing the unit safely. More workers may be involved but we expressly advise against using fewer.

 **NOTES**  
In case of doubts and/or any need for clarification during installation, contact the technical service department by e-mail, using the address: [assistenza@chinesport.it](mailto:assistenza@chinesport.it) and indicating the page and elements for which you need further details - if necessary attach a photo of the details.

**Chinesport is NOT responsible for anchoring to the walls, which must be done in a workmanlike manner, in compliance with current regulations, and entirely for the purchaser's account**

## 2 ASSEMBLY AND INSTALLATION

 **WARNING!** 

Correct installation and compliance with the loadbearing capacities are fundamental for the safety of the patient and the professional.

 **WARNINGS!** 

Before proceeding check the suitability of the environment and the structural adequacy of the available space

- The manufacturer cannot accept any responsibility for the installation done outside of their control.
- Technical skills and suitable equipment are required for a proper installation
- Unpack and handle the items with great care keeping them divided by type: if you notice any anomalies in the packaging or in the product do not continue the installation but promptly notify Chinesport SpA
- During the assembly use caution and always wear personal protective equipment
- Refer always to the installation schemes
- The manufacturer can not be responsible for the installation carried out of his control.
- Do not leave the device in use to minors.
- The equipment and accessories are intended only for the use established in the appropriate instruction manual supplied with the equipment

 **CAUTION** 

The equipment must not be used:

- Without making the proper assembly and testing.
- Tampering with and / or changes on the original product.
- In improper configuration or when there is suspicion that may cause harm to the patient.
- By children and minors unattended.
- Under weathered.

In particular:

- The installation must be performed by a competent person who has consulted the manual and found the suitability of the premises.
- Assessing hazards resulting from the association of other equipment to the structure

## 2.1 LIST OF PARTS

<b>04550</b>	<p>[A] WALL BARS          [B] FOLDING MAT          [C] BACKREST BOARD          [D] "BACKPAIN MANUAL" AND DVD          [E] MOUNTING HARDWARE</p>


## 3 ASSEMBLY

- 1 Remove the equipment from the packaging
- 2 Open the hardware bag [E] and check that all the hardware match with the quantities and types as shown in the following picture:

	[E1]	Bracket	4 pcs
	[E2]	Wood screw	8 pcs
	Spacers 6 pcs [not shown in the pictures]		

If the hardware doesn't match with the picture please contact our assistance service department

- 3 Tighten the brackets to the wallbar using the wood screws.
- 4 If the wall has a skirting, place the 6 spacers in the holes on the back of the wallbars
- 5 Check the type of the wall and its suitability to fix the wallbars.
- 6 Choose the related dowel and fix the wallbars to the wall
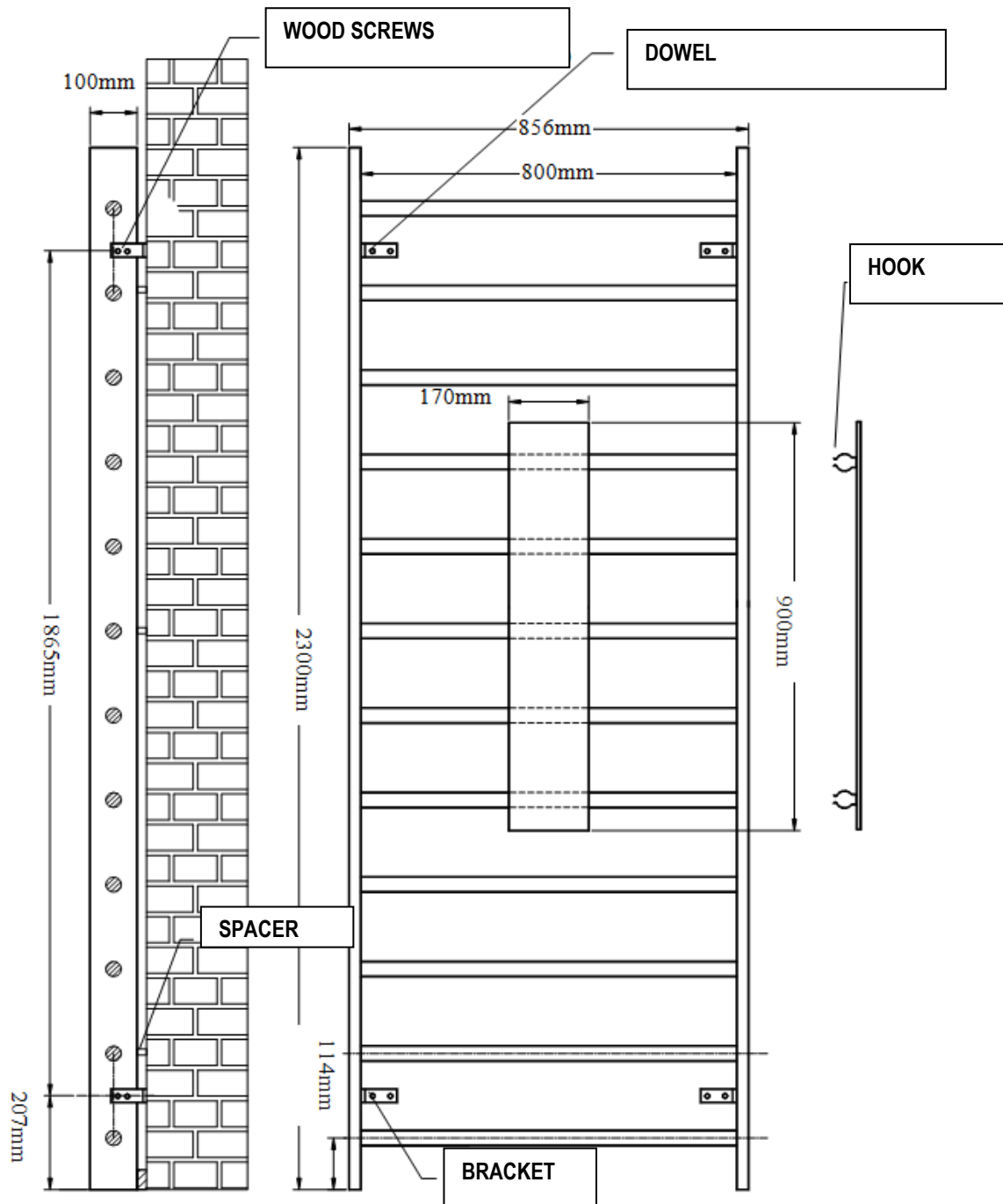
## 4 MAINTENANCE

Only qualified personnel may carry out maintenance on the equipment. Technicians must be familiar with the model, use and maintenance of the equipment.

Check the equipment whenever it has undergone a "rough treatment" (e.g.: fall, liquid spill, etc.) and whenever any doubts arise concerning its safety.

### 4.1 ROUTINE MAINTENANCE

Routine maintenance tasks are to be carried out at least once every 12 months.

FREQUENCY	TASK
Half-yearly	General cleaning of the equipment
	Visual and functional inspection of the equipment.

---

Mechanical strength, tightening the screws.
---

## **5 CLEANING AND DISINFECTION**

### **5.1 WASHING**

Dust gently and in case of light stains pass the surface with a damp cloth or a neutral detergent, rinse with water without squeezing. The following instructions are recommended for small located and more obstinate stains (oil, fat, cosmetics, ink, coffee, liquor, chewing-gum etc.): dissolve immediately with clear ethylic alcohol diluted in water to 20%, dab and then treat with a neutral detergent solution. Rinse abundantly.

### **5.2 DISINFECTION**

For the disinfection of surfaces such as chassis, sky, hand control, etc. we recommend the use of disinfectants such as AMUCHINA® 10% or equivalent disinfectant with sodium hypochlorite 0.1%

## **6 WARRANTY**

### **GENERAL AFTER-SALES SERVICE CONDITIONS FOR CHINESPORT ITEMS**

The warranty shall be valid for a period of 24 months from the date of delivery. The warranty does not cover consumable items, rechargeable batteries, and generally all the material subject to wear, breakdowns caused by knocks, falls, misuse or improper use, accidental events and transport damage. This warranty shall be automatically invalidated if the equipment is tampered with. Repair of the equipment will be carried out in relation to the breakdown stated. An estimate must be specifically requested. Payment is cash on delivery, unless other agreements are made. Any dispute will be dealt with solely by the Court of UDINE.

### **WARRANTY REPAIRS**

Repairs under warranty must be specifically requested. Repairs under warranty shall be carried out at our workshop and upon authorization, complete with the return authorization number. The return shipment for products sent in their original package will be made free of carriage charges. A fiscal document (bill of sale, purchase invoice or cash receipt) is required to make a warranty claim; the purchase date indicated must fall within the warranty time period. If the product is found to be in good working order without defect, the fault-finding service shall be charged to the customer.

### **OUT-OF-WARRANTY REPAIRS**

For out-of-warranty repairs carried out at our workshops upon authorization, transportation charges shall be paid by the customer. The warranty of the parts replaced during the repair operations is 12 months from the moment the equipment is collected. If the product is found to be in good working order without defect, the fault-finding service shall be charged to the customer.

### **REPAIRS MADE AT HOME**

In the event repairs are made at the customer's home, a written request is required indicating the complete details of the person making the request, the type of machine and the breakdown. The cost per kilometer for on-site technician visits will be agreed upon in relation to the urgency of the customer. In the event the machine is under warranty, the customer will have to pay only the expenses incurred for the technician to reach the customer's house. Time is calculated from the moment the technician leaves our workshop until he/she returns. The time of his/her return will be estimated based on the time required to get there.

### **SALE OF SPARE PARTS AND CONSUMABLE ITEMS**

The orders must be in writing and specify the customer, the item and the shipping method as well as the precise fiscal data of the customer. The order processing time can vary depending on the availability of the material. Returned spare parts will not be accepted. Payment will be cash on delivery, unless other agreements are made.

## 7 MAINTENANCE RECORD SHEET

EQUIPMENT		DATE of first installation	SERIAL NUMBER		
DATE OF SERVICE	OPERATIONS CARRIED OUT	TECHNICIAN	SIGNATURE	NEXT VERIFY DATE	

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