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USER MANUAL AXILLARY WALKER

INTRODUCTION

Thank you for choosing Demarta Virginio. The products have been designed and built to meet all your needs for a practical, correct and safe use. This manual contains some small suggestions for a correct use of the device of your choice and some precious advice for your safety. We recommend that you carefully read the entirety of this manual before using the item you purchased. If in doubt, please contact your dealer, who will be able to help and advise you correctly.

The manufacturer reserves the right to make changes to the device and to this manual without prior notice in order to improve its characteristics.

GENERAL WARNINGS

- Use the product in the manner indicated in this manual. If doubts persist, contact DEMARTA VIRGINIO S.r.l..
- Always evaluate the compatibility of the product with the characteristics of the subject together with the doctor or therapist.
- Always check the integrity of the product before each use, even when new. Do not use if damaged or worn.
- DEMARTA VIRGINIO S.r.l. declines any responsibility for damage to property or persons resulting from improper use, not agreed with a professional figure, not corresponding to what is indicated in the user manual.
- Any accidents, anomalies, breakages must be reported to DEMARTA VIRGINIO S.r.l. in order to monitor the quality and criticalities of the product.

DECLARATION OF CONFORMITY

DEMARTA VIRGINIO S.r.l. declares under its sole responsibility that the product made and traded by DEMARTA VIRGINIO S.r.l. comply with the applicable provisions of the 2017/745 Regulation on Medical Devices of 5 April 2017.

For this purpose, DEMARTA VIRGINIO S.r.l. guarantees and declares under its sole responsibility what follows:

- The devices satisfy the requirements of general safety and performance requested by the Annex I of regulation 2017/745 as laid down by the Annex IV of the above mentioned regulation.
- The devices ARE NOT MEASURING INSTRUMENTS.
- The devices ARE NOT MADE FOR CLINICAL TESTS.
- The devices are packed in NON-STERILE BOX.
- The devices belong to class I in accordance with the provisions of Annex VIII of the above mentioned regulation.

INTENDED USE

Walkers are indicated in the phase of re-education to walking in elderly subjects or with functional limitations of various kinds.

They offer stability, ease of use and improve the user's sense of balance.

This aid is designed for use by a single user in indoor environments.

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- The walker should only be used by patients with adequate motor skills that ensure its safe use.
- Consult your doctor or therapist to determine the suitability of the product, the proper adjustment and the correct use of the walker.
- Your doctor or therapist should assist in the height adjustment of the walker.
- The walker can only be used on a firm and flat surface.
- Do not use the walker to climb stairs.
- Pay particular attention to wet surfaces.
- All wheels must be in contact with the floor at all times during use to ensure proper balance of the walker and maximum safety.

Before using the device always check:

- That all parts of the walker are securely locked.
- Make sure the seat is positioned correctly, if present.
- Check that the brakes are working correctly.

DESCRIPTION

The Anti-Brachial Walker has a foldable and height-adjustable steel structure. The height adjustment occurs in 6 points, ensuring excellent customization to adapt to various heights and preferences. The two underarm supports, also adjustable, ensure stable and comfortable support during use.

The walker is equipped with four 10 cm swivel wheels, which allow great maneuverability, making it particularly suitable for moving even in narrow spaces. The rear wheels are equipped with a brake, which offers greater safety and allows you to block the walker when necessary, ensuring stability during breaks or when you are stationary.

TECHNICAL DATA

Code	Length	Inside Width	Width	Height	Seat Dimension	Floor to Seat Height	Floor to Support Bar Height	Weight Capacity	Wheel Dimension	Net Weight	Carton Size
Measures in Cm											
01.06.650S	74	46	51	120/165	43x21	44.5/57	99/113	120	10 (4")	13	73x11x80

ASSEMBLY

1) Remove the frame from the packaging and open it by pushing the two side supports outwards.

2) Remove the frame connecting U-tube from the box.

• Insert the U-shaped tube into the appropriate holes at the two ends of the frame

- Press the metal spring button so that it retracts completely into the tube.

- Release the metal spring button in front of the exact hole position.

⚠ ATTENTION!

Repeat the operation performed on the other end of the frame.

3) Remove the 4 wheels from the packaging.

⚠ ATTENTION!

Insert the two wheels with brakes into the rear of the frame and the two wheels with brakes into the front of the frame.

- Press the metal spring button so that it retracts completely into the tube.
- Select the desired height by adjusting the legs with an upward or downward movement.
- Release the metal spring button in front of the desired hole position.

⚠ ATTENTION!

Before use, always make sure that each spring button protrudes fully from the selected hole position and that the height is the same for both 4 wheels.


- Always keep the braked wheels locked during assembly by pushing the grey pedal with your foot.
- To unlock the wheels, return the brake upwards to the horizontal position.

4) Remove the pair of armpit supports from the packaging.

Insert the two supports into the two holes provided at the two ends of the frame.

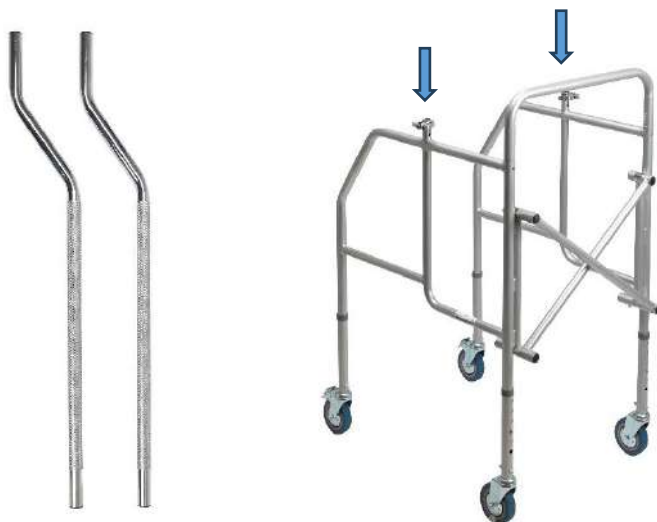You can adjust the height of the armpit supports with an upward or downward movement. After choosing the desired height, close the locking lever in the knurled area on the armpit supports.

Insert the two armpit supports into the two holes provided at the two ends indicated.

OPTIONAL

It is possible to insert the **Deluxe Saddle** onto the walker frame.

Place it over the frame and apply light pressure on it until the 4 clips lock correctly.

USE

- Start pushing when you are standing with both hands on the handles.
- Make sure your hands are dry before using the walker to avoid the risk of slipping.
- Walk upright and look straight ahead, not down.
- Do not use the walker to negotiate steps.

CLEANING

To clean the parts, use only a damp cloth with neutral detergents. Do not use detergents containing solvents and metal brushes or abrasive or sharp objects, they could damage the product. Absolutely avoid the use of acidic, alkaline products or solvents such as acetone or thinner.

GENERAL DISPOSAL CONDITIONS

When disposing of the device, never use housing solid waste collection systems.

It is necessary to dispose of the device through the common ecological islands common for the recycling of used materials.

WARRANTY

DEMARTA products are guaranteed from material or manufacturing faults for 2 years from the purchasing date.

DEMARTA does not guarantee for damages caused by the following condition: natural disaster, not authorized maintenance and failure to adhere to the guidelines and instructions for use.